

File D:\DATA\DCSG\MONITOR\P2-00418-LCMSA048.D

Injection Date : 7 Mar 12 10:52 am +0800

Tgt Mass(EZX) :

Sample Name : P2

Location : P1-B-02

Acq. Operator : DCSG-monitor

Inj : 1

Spec. Reported : MS Integration

Inj Volume : 2 ul

Acq. Method : D:\METHODS\ROUTINE.M

Analysis Method : D:\METHODS\ROUTINE.M

Sample Info : *** No target masses specified ***

Method Info : Column:Waters Sunfire C18(30mm*4.6mm,2.5um)

Detection: UV at 214 nm, Flow rate: 1.5 mL/min


Eluent A1: 0.025 % TFA in Water, Eluent B1: 0.025 % TFA in ACN

Run Time: 5.0 min


Gradient: Time (min) 0 3.0 5.0

B1 (%) 10 95 95


ELS1 A, Voltage

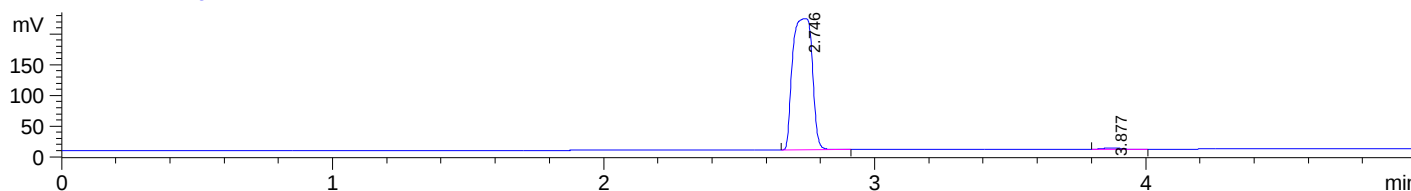


DAD1 A, Sig=214,4 Ref=off

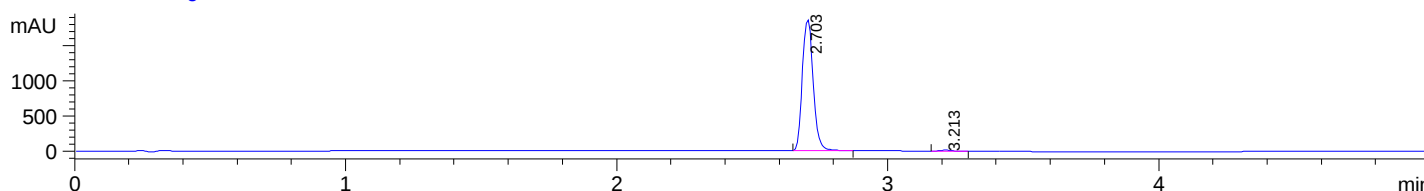

DAD1 B, Sig=254,16 Ref=off

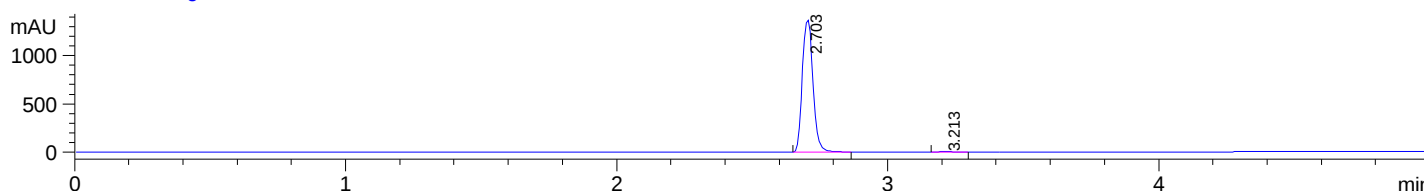

MSD1 TIC, MS File ES-API, Pos, Scan, Frag: 70

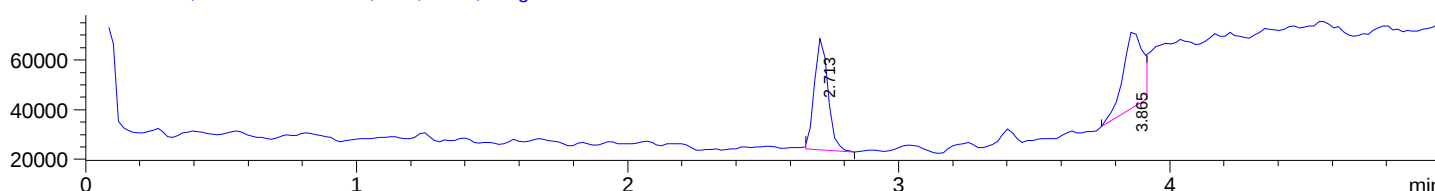

Integration Results for ELS1 A, Voltage

RetTim	Width	Area	Height	Area%
2.75	0.08	1072.78	213.67	99.33
3.88	0.06	7.22	1.43	0.67

Integration Results for DAD1 A, Sig=214,4 Ref=off

RetTim	Width	Area	Height	Area%
2.70	0.04	5251.03	1867.80	99.38
3.21	0.04	32.80	12.41	0.62

Integration Results for DAD1 B, Sig=254,16 Ref=off


RetTim	Width	Area	Height	Area%
2.70	0.04	3799.14	1376.53	99.44
3.21	0.04	21.39	7.98	0.56

Integration Results for MSD1 TIC, MS File

2.71	0.06	157106.61	44965.88	50.70
3.87	0.07	152739.38	30973.94	49.30

RetTim	Width	Area	Height	Area%
--------	-------	------	--------	-------

Ret. Time: 2.71 <<<< POSITIVE SPECTRA >>>>

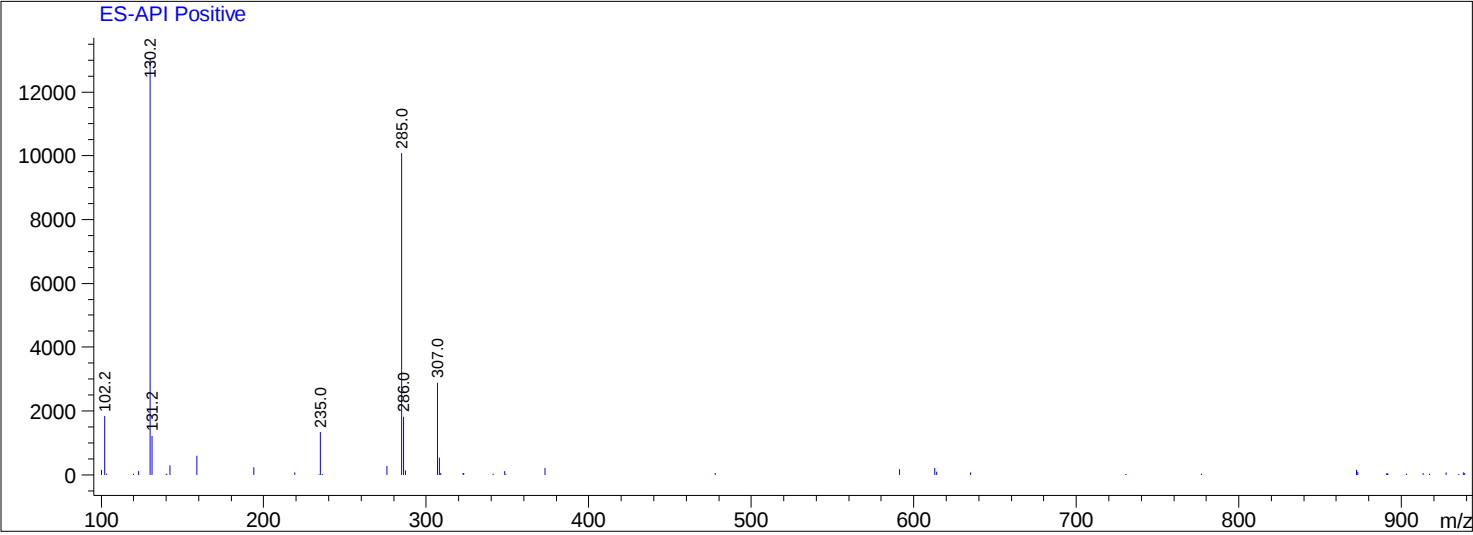

Ret. Time: 3.87

