


KI0673 HPLC Analysis Report

Sample: Cefepime
 Sample ID: KI0673
 Buffer A: 0.1% TFA in 100% water (v/v)
 Buffer B: 0.1% TFA in 80% acetonitrile + 20% water (v/v)
 Gradient: 10-70% Buffer B in 20min
 Flow: 1ml/min Wavelength: 220nm
 Column: Venusil C18 5um 4.6*250mm

Chromatogram

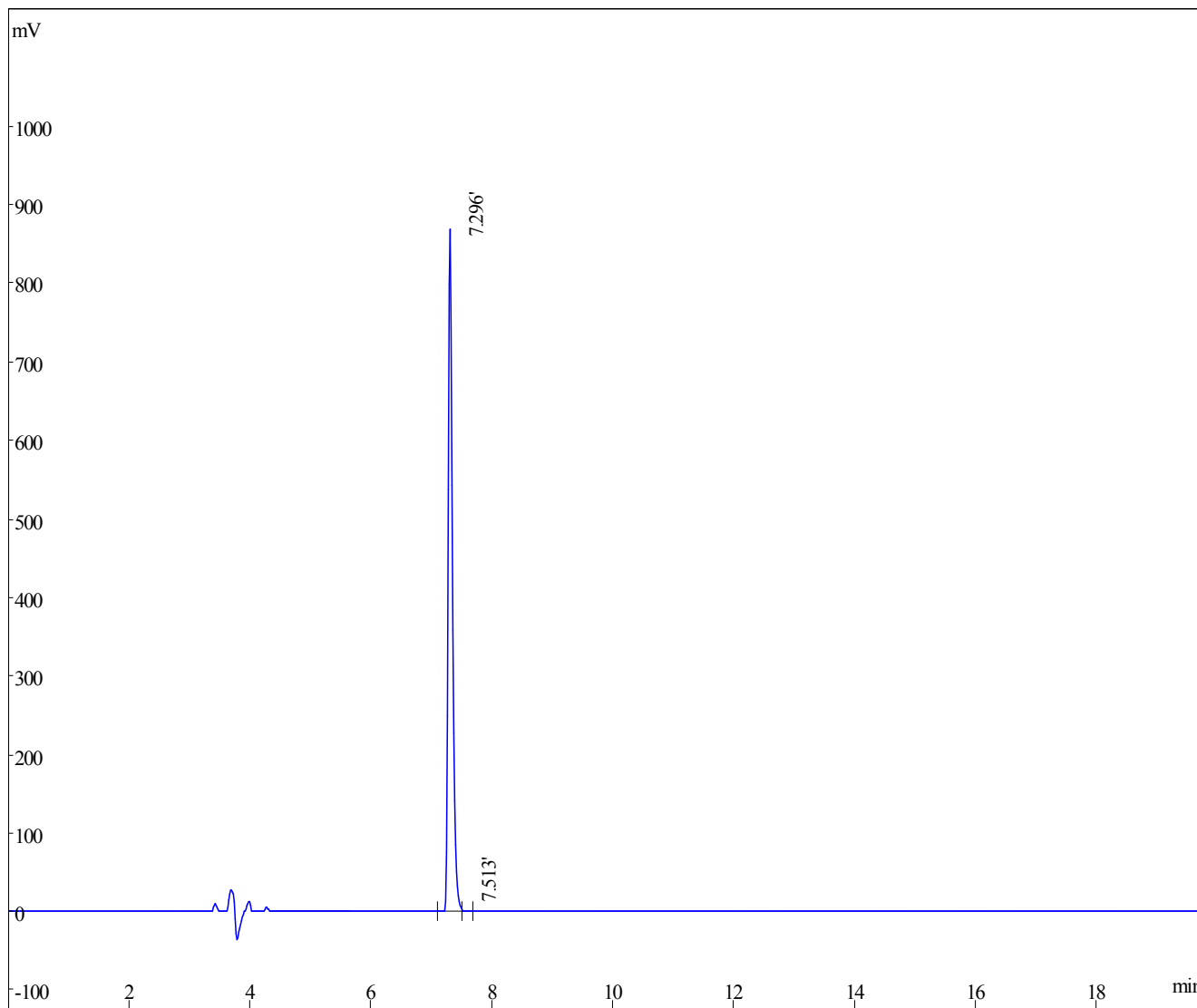

Signal	Peak#	RT (min)	Area%	Area (mAU*s)	Height (mAU)
	1	7.296	99.78	4301259	882125
	2	7.513	0.222	9571	1944
Total			100	4310830	884069