

KI4446 HPLC Analysis Report

Sample: DMAB-anabaseine dihydrochloride

Sample ID: KI4446

Buffer A: 0.1% TFA in 100% water (v/v)

Buffer B: 0.1% TFA in 80% acetonitrile + 20% water (v/v)


Gradient: 10–100% Buffer B in 30min

Flow: 1ml/min

Wavelength: 220nm

Column: Venusil C18 5um 4.6*250mm

Chromatogram

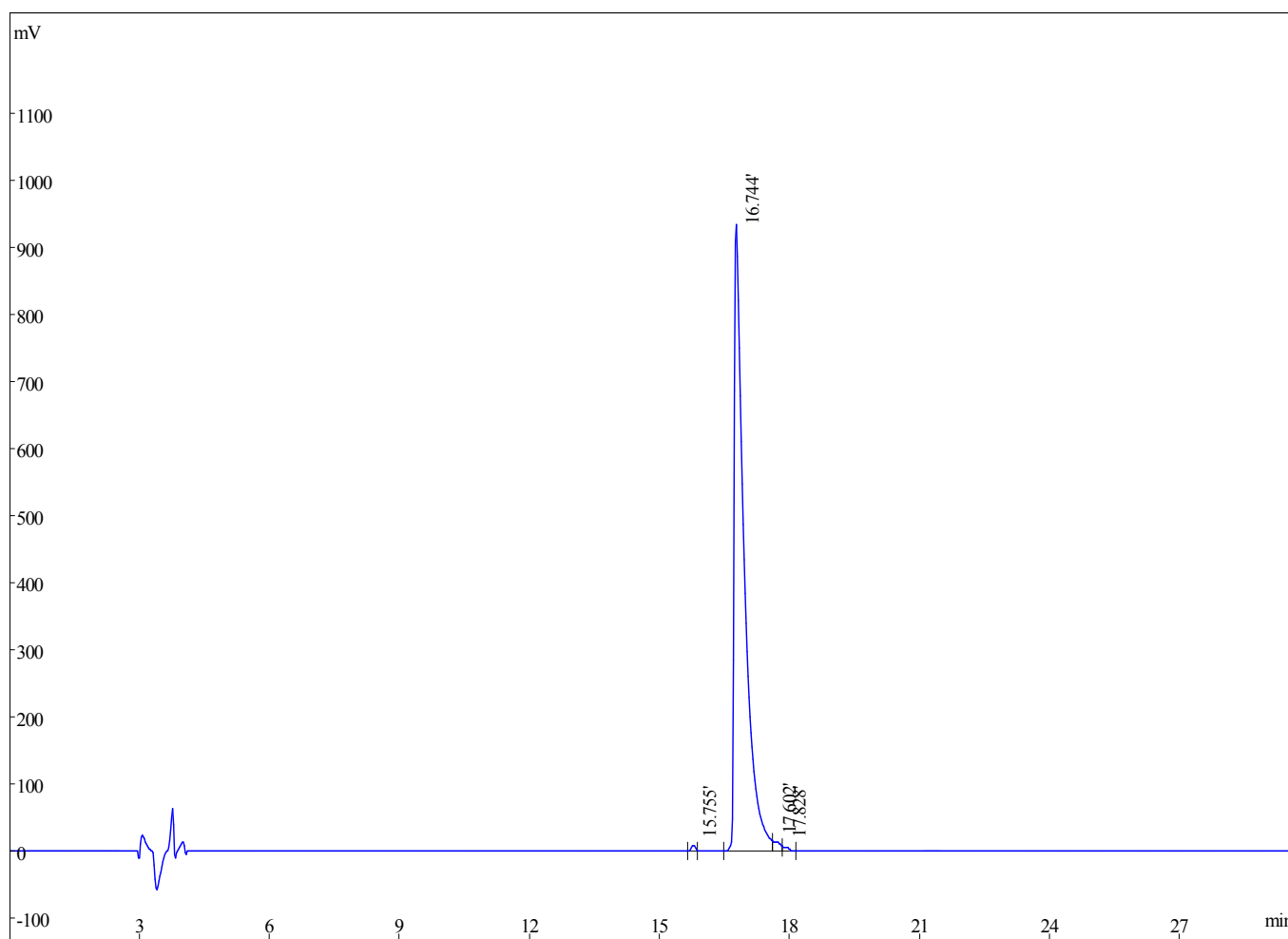

Signal	Peak#	RT(min)	Area%	Area(mAU*s)	Height (mAU)
	1	15.755	0.4303	66926	10030
	2	16.744	98.02	15246229	941595
	3	17.602	1.094	170146	14804
	4	17.828	0.4589	71376	7426
Total			100	15554677	973855