

KI0927 HPLC Analysis Report

Sample: Decitabine

Sample ID: **KI0927**

Buffer A: 0.1% TFA in 100% water (v/v)


Buffer B: 0.1% TFA in 80% acetonitrile + 20% water (v/v)

Gradient: 2% Buffer B in 10min

Flow: 1ml/min Wavelength: 220nm

Column: Venusil C18 5um 4.6*250mm

Chromatogram

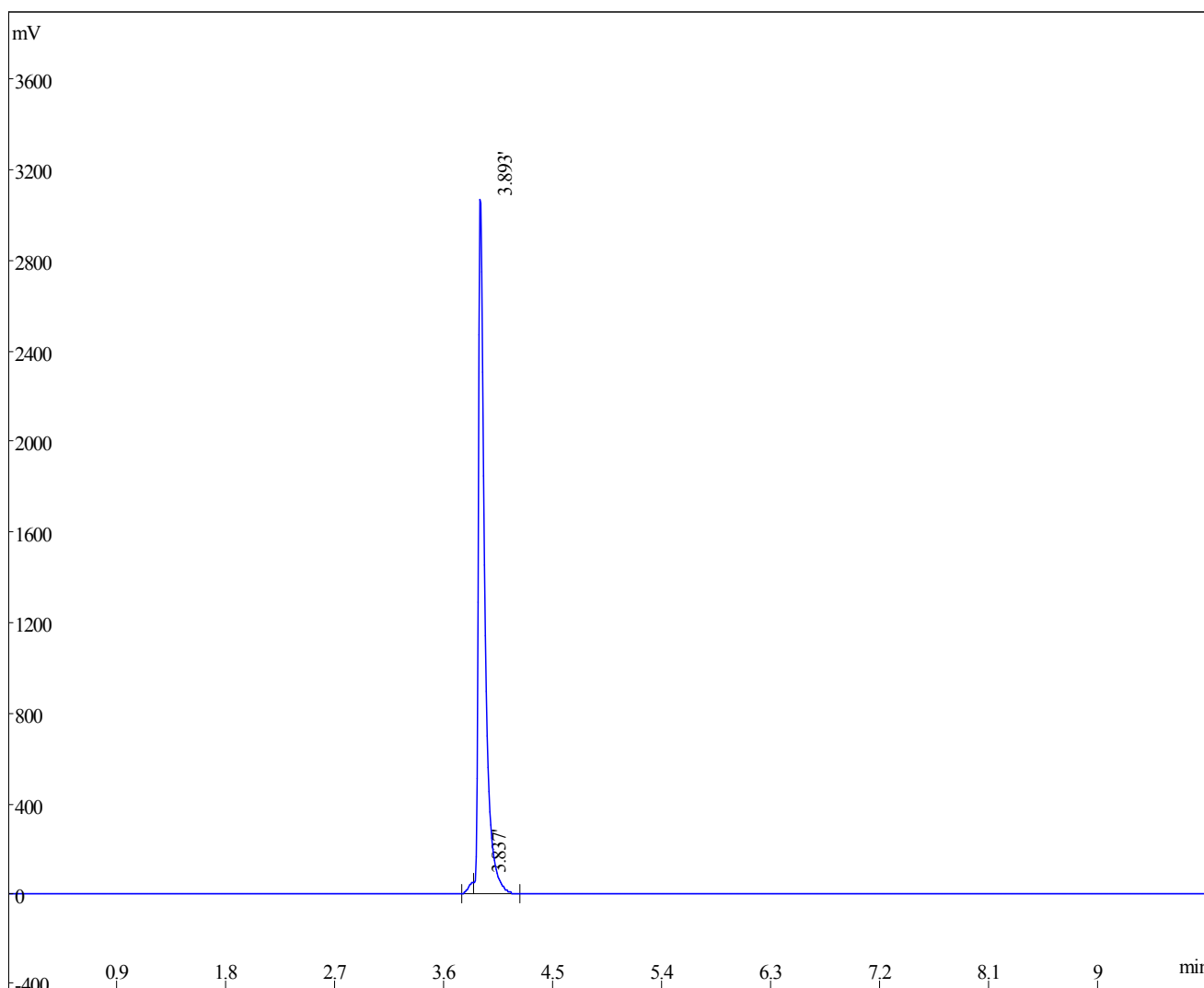

Signal	Peak#	RT (min)	Area%	Area (mAU*s)	Height (mAU)
	1	3.837	1.528	175404	52121
	2	3.893	98.47	11304597	3120424
Total			100	11480001	3172545