

HPLC Analysis Report

Sample: Finasteride

CAS No.: 98319-26-7

Batch No.: KB20140521

Mfg. Date: May 21, 2014

Method Info: 10%-80% Buffer B 0-25 MIN

Buffer A: 0.04% TFA in 100% water(v/v)

Buffer B: 0.04% TFA in 100% acetonitrile (v/v)


Gradient: 10-80% Buffer B in 25 min

Column: Alltima C18, 4.6*250mm, 5µm

Flow Rate: 1mL/min

Wavelength: 220 nm

Date: May. 27, 2014

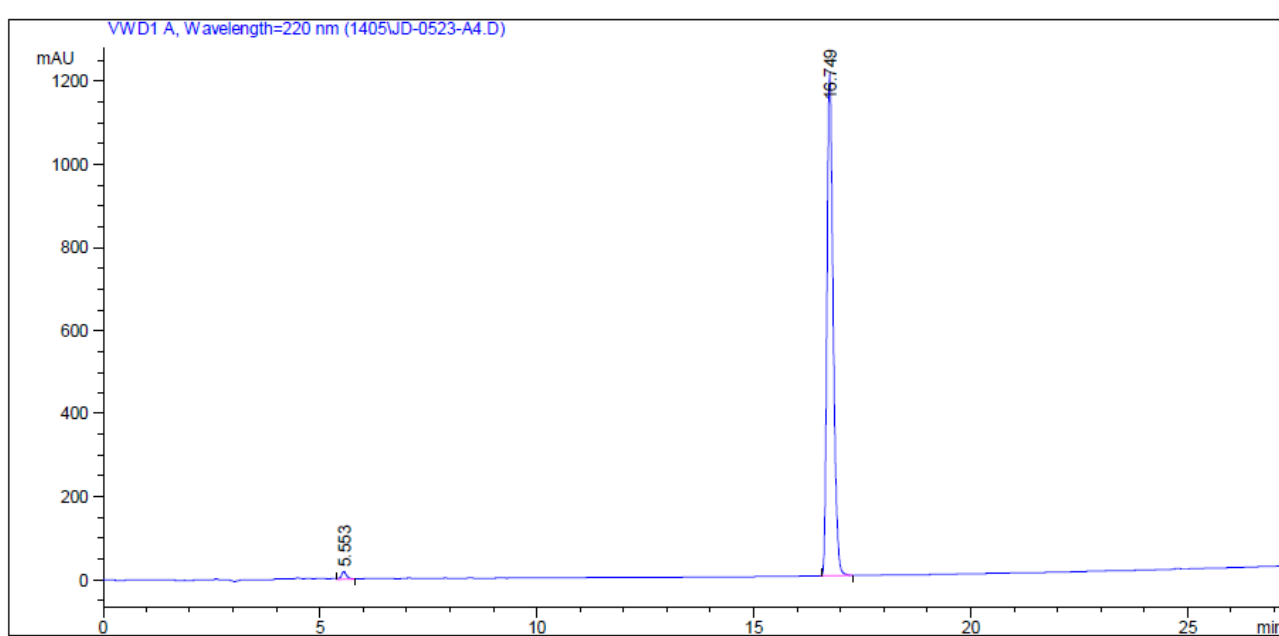

Peak#	RetTime(Min)	Width(Min)	Area(mAU*s)	Height(mAU)	Area(%)
1	5.553	0.1250	149.11472	18.13680	1.2284
2	16.749	0.1512	1.19902e4	1203.09387	98.7716
Totals			1.21393e4	1221.23067	100

Analysis Result:

Purity: 98.77%