

KI0483 HPLC Analysis Report

Sample:LDE225

Sample ID: **KI0483**

Buffer A:0.1% TFA in 100% water(v/v)

Buffer B:0.1% TFA in 80% acetonitrile + 20% water(v/v)


Gradient:65-85% Buffer B in 20min

Flow:1ml/min

Wavelength:220nm

Column:Venusil C18 5um4.6*250mm

Chromatogram

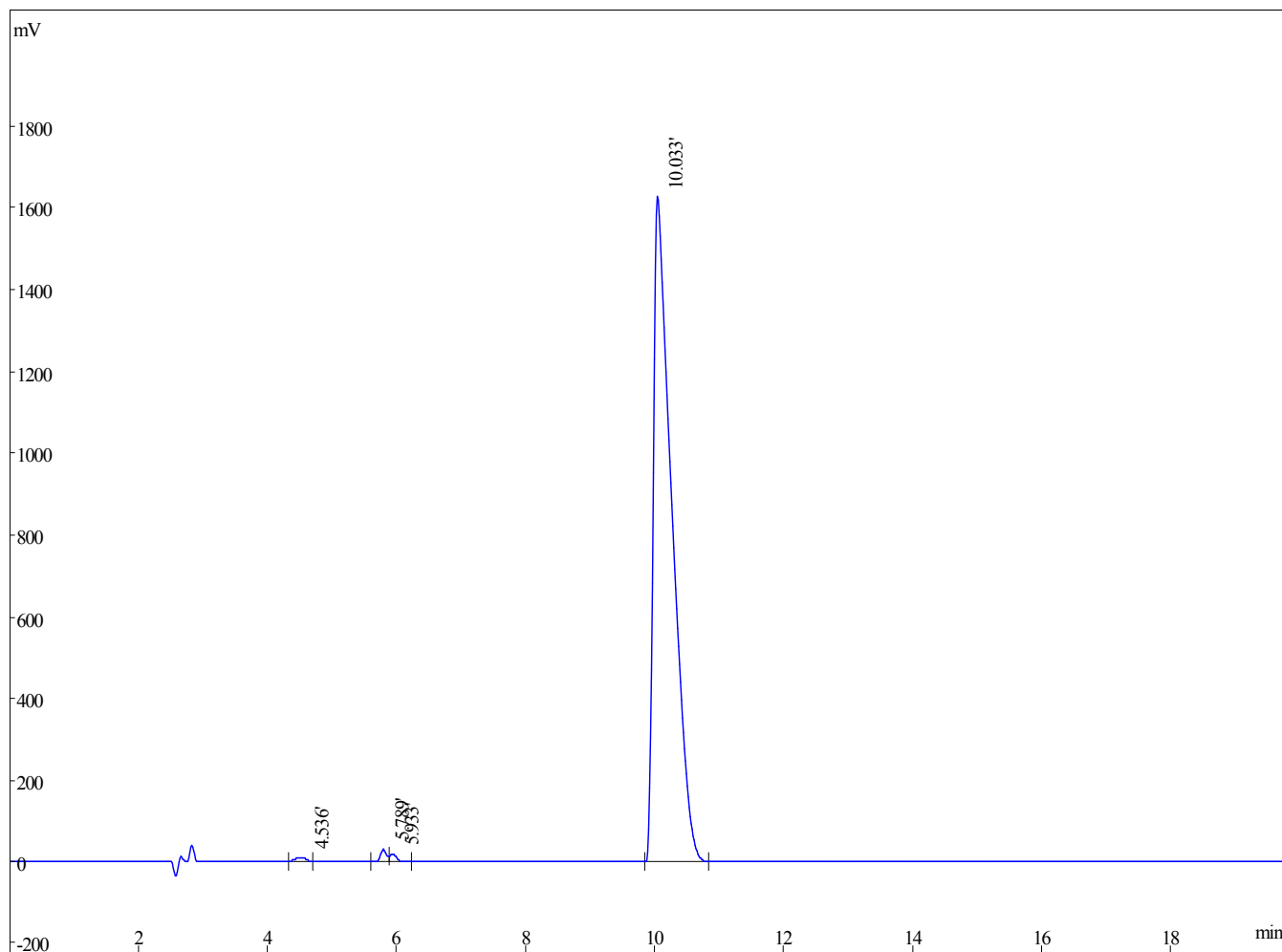

Signal	Peak#	RT(min)	Area%	Area(mAU*s)	Height (mAU)
	1	4.536	0.4717	154707	12801
	2	5.789	0.7042	230962	30905
	3	5.933	0.4029	132135	19787
	4	10.033	98.42	32280128	1631672
Total			100	32797932	1695165