

HPLC Analysis Report

Sample: Cariporide

Lot No.: KB20150527

Buffer A: 0.04% TFA in 100% water(v/v)

Buffer B: 0.04% TFA in 100% acetonitrile (v/v)

Gradient: 10-80% Buffer B in 25 min then 100% Buffer B for 5 min.


Flow: 1ml/min

Wavelength:220nm

Column: Alltima C18,4.6*250mm,5um

Date: May 27, 2015

Chromatogram

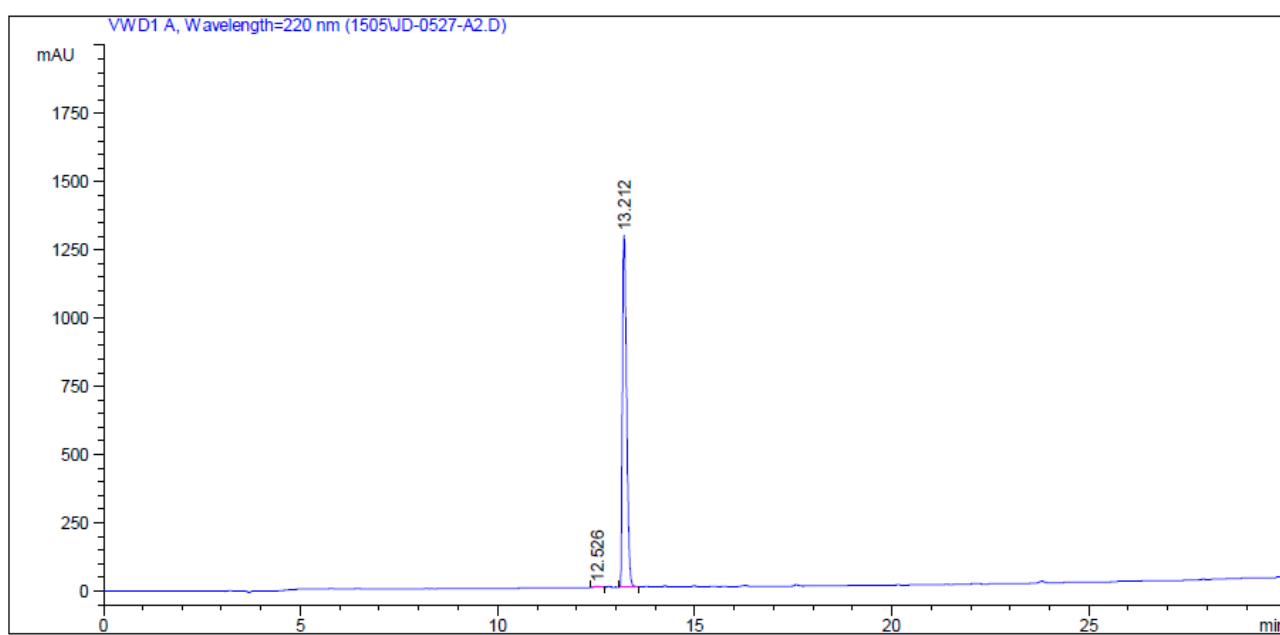

Peak#	RetTime(Min)	Width(Min)	Area(mAU*s)	Height(mAU)	Area(%)
1	12.526	0.1665	60.10826	5.33119	0.6321
2	13.212	0.1147	9448.83008	1288.32507	99.3679
			9508.93834	1293.65626	100

Analysis Result:

Purity: 99.37%