

MK-677 HPLC Analysis Report

Sample: **MK-677**

Sample ID: **KI3636**

Buffer A: 0.1% TFA in 100% water (v/v)

Buffer B: 0.1% TFA in 80% acetonitrile + 20% water (v/v)


Gradient: 10-100% Buffer B in 30min

Flow: 1ml/min

Wavelength: 220nm

Column: Venusil C18 5um 4.6*250mm

Chromatogram

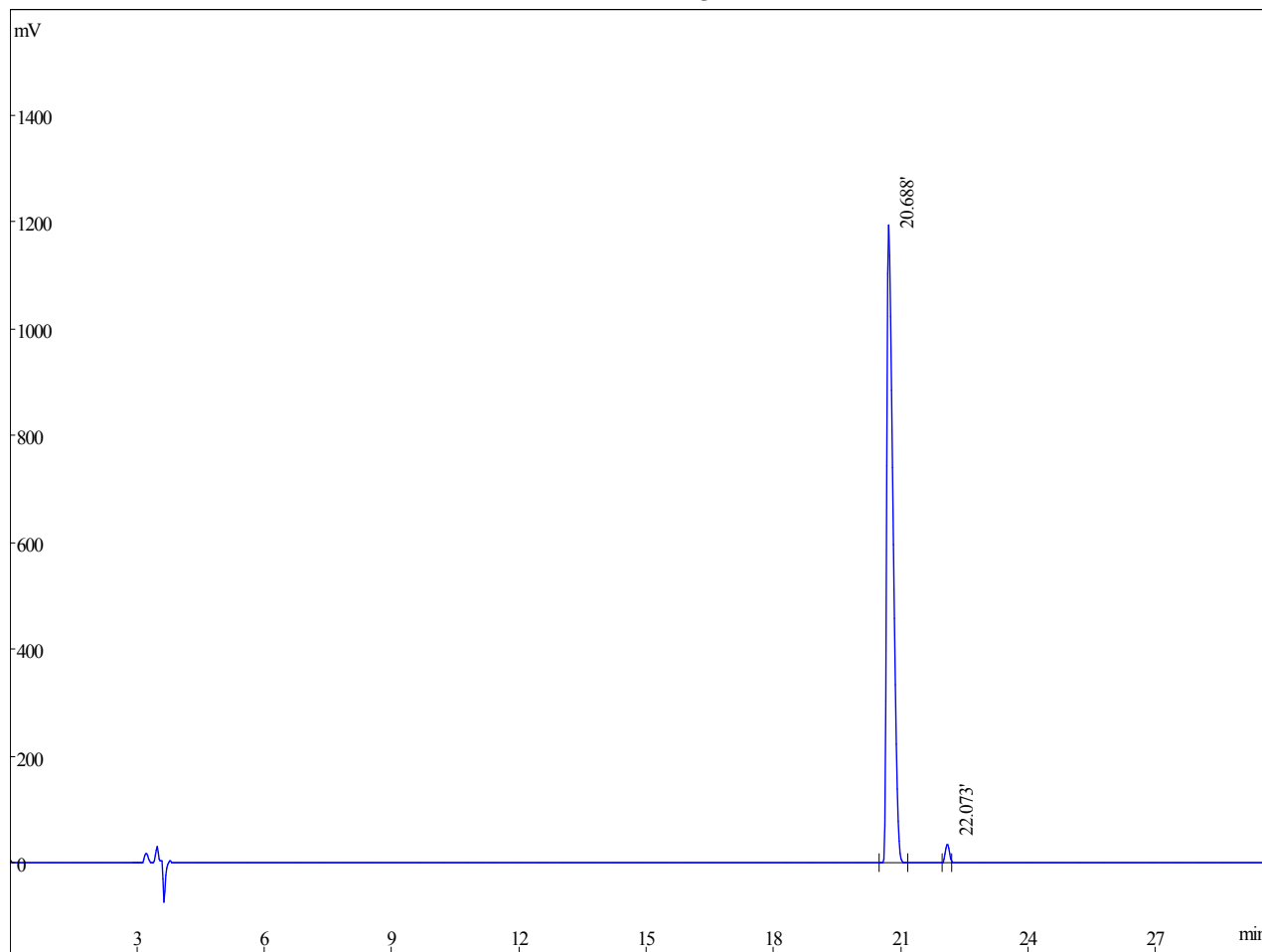

Signal	Peak#	RT(min)	Area%	Area(mAU*s)	Height (mAU)
	1	20.688	98.05	12316278	1199751
	2	22.073	1.949	244752	35775
Total			100	12561030	1235526