

## KI0222 HPLC Analysis Report

Sample: MLN8237

Sample ID: KI0222

Buffer A: 0.1% TFA in 100% water (v/v)

Buffer B: 0.1% TFA in 80% acetonitrile + 20% water (v/v)


Gradient: 59–79% Buffer B in 20min

Flow: 1ml/min

Wavelength: 220nm

Column: Venusil C18 5um 4.6\*250mm

### Chromatogram

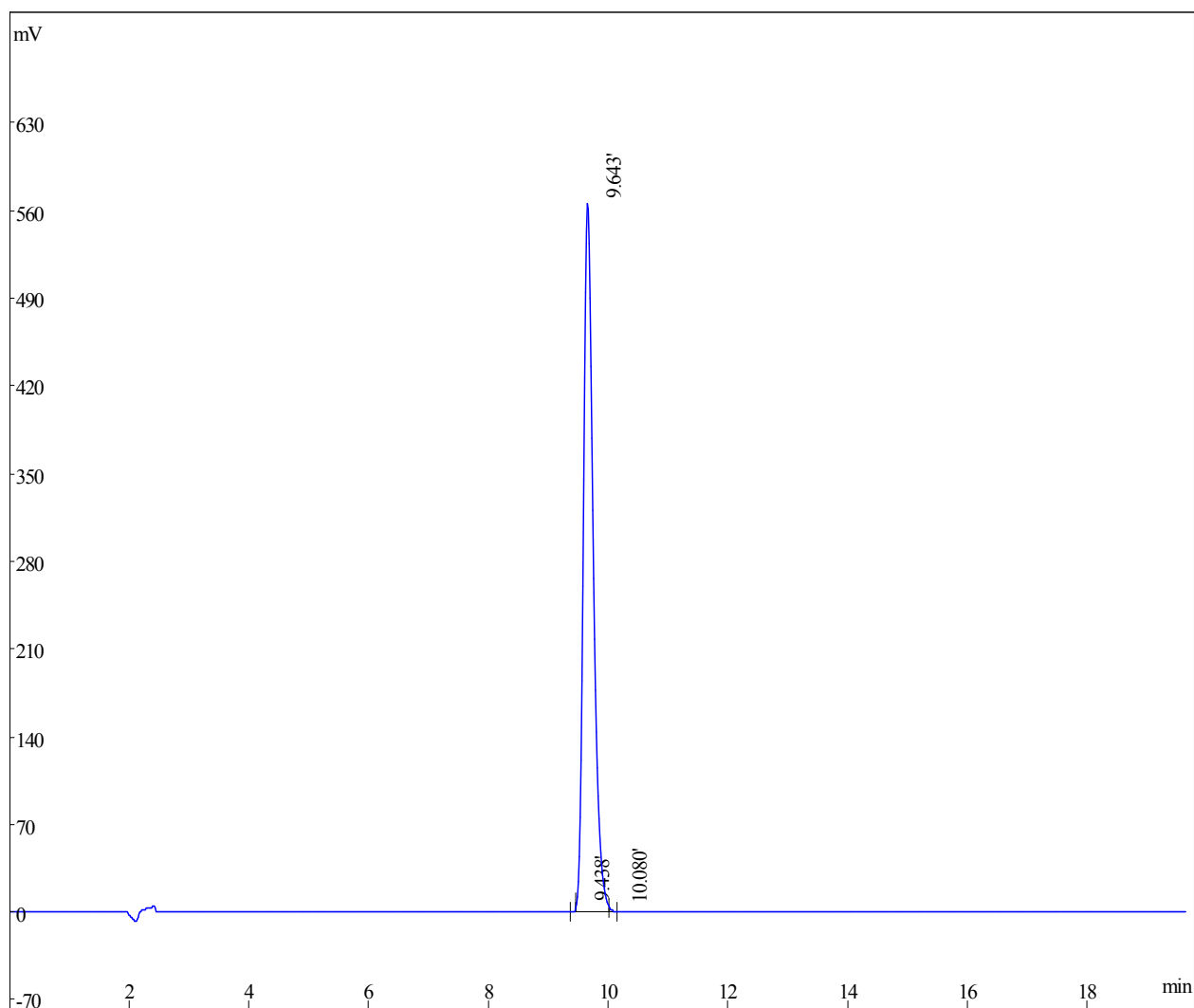

| Signal | Peak# | RT(min) | Area%  | Area(mAU*s) | Height (mAU) |
|--------|-------|---------|--------|-------------|--------------|
|        | 1     | 9.438   | 0.1448 | 9380        | 2276         |
|        | 2     | 9.643   | 99.68  | 6456068     | 567636       |
|        | 3     | 10.080  | 0.1703 | 11026       | 1076         |
| Total  |       |         | 100    | 6476474     | 570988       |